

PRODUCT INFORMATION




Goat Anti-Human IgG (H+L) HRP

Item No. 46750

Overview and Properties

Contents:	This vial contains 100 µg of immunoaffinity-purified IgG conjugated to HRP.
Synonyms:	Horseradish Peroxidase-Conjugated Goat Anti-Human IgG (H+L)
Cross Reactivity:	(+) IgG
Species Reactivity:	(+) Human
Form:	Liquid
Storage:	-20°C (as supplied)
Stability:	≥1 year
Storage Buffer:	15 mg/ml BSA, pH 7.6, with 0.01 M sodium phosphate, 0.25 M sodium chloride, and 50% glycerol
Host:	Goat
Applications:	ELISA and Western blot (WB); the recommended starting dilution is 1:10,000. Other applications were not tested, therefore optimal working concentration/dilution should be determined empirically.

Image

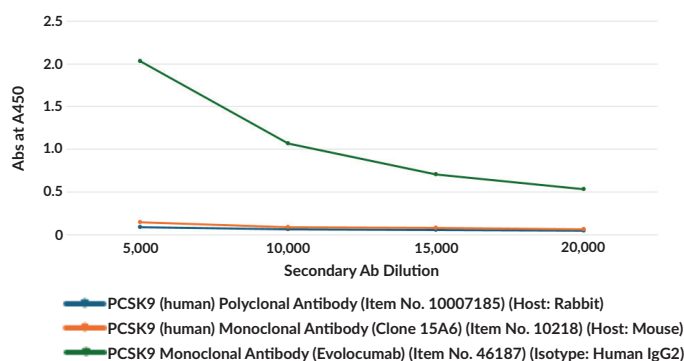

ELISA. Plates were coated with PCSK9 (human, recombinant) (Item No. 20631), blocked with BSA, and probed with the indicated primary antibody at the recommended dilutions and Goat Anti-Human IgG (H+L) HRP (Item No. 46750) secondary antibody with TMB substrate and an H_2SO_4 stop solution.

Description

Goat Anti-Human IgG (H+L) HRP is a polyclonal secondary antibody that can be used to detect human IgG in ELISA and Western blot (WB) applications. This HRP-conjugated secondary antibody can be quantified using various substrates, including TMB (Item Nos. 36361 | 70450).

WARNING
THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA
This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the complete Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY
Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 10/23/2026

CAYMAN CHEMICAL
1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA
PHONE: [800] 364-9897
[734] 971-3335
FAX: [734] 971-3640
CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM